



# Lavieen

Lavieen (Thulium Laser/1,927nm), well known as BB Laser, shows significant brightening outcome. It is a non-ablative skin resurfacing and rejuvenation laser delivering remarkable results with minimal downtime.

## Lavieen

Thulium Laser / 1,927nm

Irradiation Type	Thulium
Wavelength	1,927nm
Diameter	300 $\mu$ m ( 10, 20, 30mm tips)
Output	Max. 10W(Adjusting by Watt & Pulse duration)
Scan Size	1 x 1mm ~ 30 x 30mm (1 x 1mm step)
Scan Distance	0.2~2.0mm (0.1mm step)
Tips	10mm, 20mm, 30mm
Dimension	279 x 686 x 849 (W x D x H)mm
Weight	30kg
Cooling System	Air cooling
Electrical	220VAC, 50/60Hz, 800VA

**BELIUM**  
MEDICAL

Avda. Dionisio Cifuentes, 34  
33203 Gijón · Asturias · España  
**985 30 12 34**  
info@belium.es  
[www.beliummedical.es](http://www.beliummedical.es)

A9

For Export Only

## Lavieen

Thulium Laser / 1,927nm



**BELIUM**  
MEDICAL

# Get Brightened without B.B cream




Lavieen is excellent way to treat melasma, pigment lesions, wrinkles, large pore and more by adjustable energy intensity and pulse duration simultaneously during procedure.



## B.B Laser

### *Sub-ablative*

- Target to the junction between epidermis and dermis
- Coagulation and denaturalization
- Mainly for whitening skin
- Less downtime

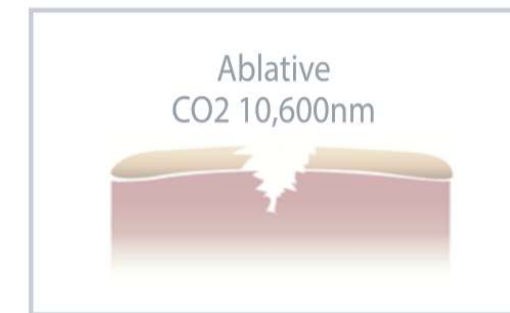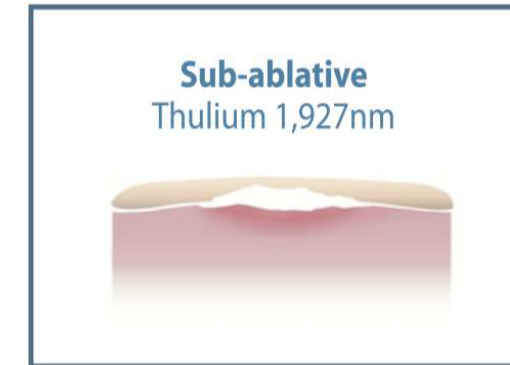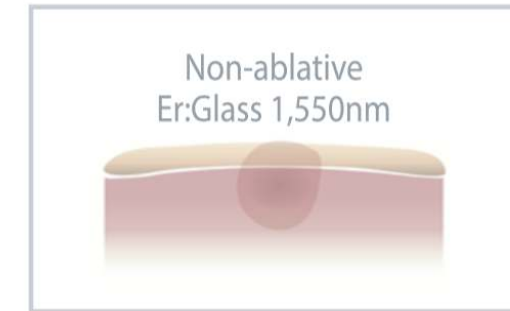

## Treat Diverse Indications

by separate setting of "Pulse duration" and "Fluence"

Low maintenance cost owing to non-consumable and long life time of thulium fiber

## Skin Whitening

- Pigmentary treatment
  - Skin tone (whitening)
  - Epidermal melasma
  - Freckles
  - Aged spots
- Dermal remodeling
  - Skin texture
  - Wrinkles
  - Lifting & tightening
  - Large pores

