



Specification

Model	Intima RF (SLM-7010E)
Performance	Maximum RF Output (150W) / Frequency (480KHz) Maximum Operation Time (Shot type: 0.5~2.5sec / Time: 15min) Impedance Check Range (1~100W)
Display	10.2 inch ColorTouch LCD
Input Power	AC100~240V, 50/60Hz
Power consumption	200W
Dimension	W34 x D458 x H1284
Weight	23kg

BELIUM
MEDICAL

DISTRIBUIDOR OFICIAL
Avda. Dionisio Cifuentes, 34
33203 Gijón 985 30 12 34
info@belium.es
www.beliummedical.es

Intima RF®

Temperature and Impedance Dual-Controlled
Radiofrequency Vaginal Rejuvenation

- Vaginal Tightening
- Labia Remodeling
- Stress Urinary Incontinence
- Skin Rejuvenation


Maximum Safety and Effectiveness through
Smart OED (Optimized Energy Delivery)
System



BELIUM
MEDICAL

Intima RF®

Temperature and Impedance Dual-Controlled Radiofrequency Vaginal Rejuvenation



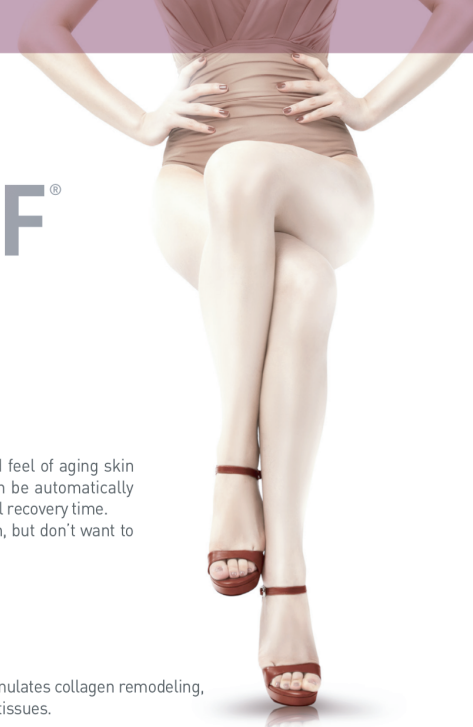


Overview

Intima RF is a non-invasive radiofrequency therapy that can address the look and feel of aging skin on many areas of body especially on the most private parts of the women. It can be automatically customized to each patient with the smart OED system, and typically require minimal recovery time. It's a perfect solution for the patients who want to tighten and rejuvenate their skin, but don't want to miss out on their work, social or sexual life.

Principle

Focused RF thermal energy evokes fibroblasts to produce new collagen fibers. It stimulates collagen remodeling, thickens the vaginal canal walls and promotes the functional recovery of the vaginal tissues.

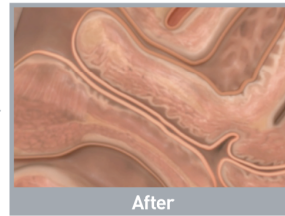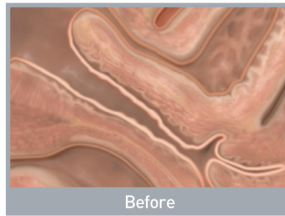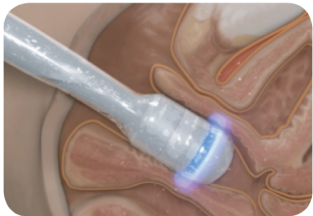

Effect



Vaginal Tightening

RF thermal effect on the vaginal mucosa, stimulates collagen remodeling, revitalizing the tissues and tightens the vaginal canal.



Labia Remodeling

Labia Remodeling using RF reduces the size and diameter of the labia or vagina lips, making them smaller, more aesthetically appealing, and asymmetrical.



Stress Urinary Incontinence

Enhanced thermal and shrinkage effect on the endopelvic fascia, ligament and mucosal tissue of the anterior vaginal canal to improve stress urinary incontinence



Skin Rejuvenation

RF energy rejuvenates facial skin making much younger, sleek and appealing look



Character

• Smart OED (Optimized Energy Delivery) System

Real-time monitoring of the skin impedance and temperature of the each patient, enabling highly effective and safe procedure.

• Versatile Function

Multiple handpiece options enabling multiple treatments

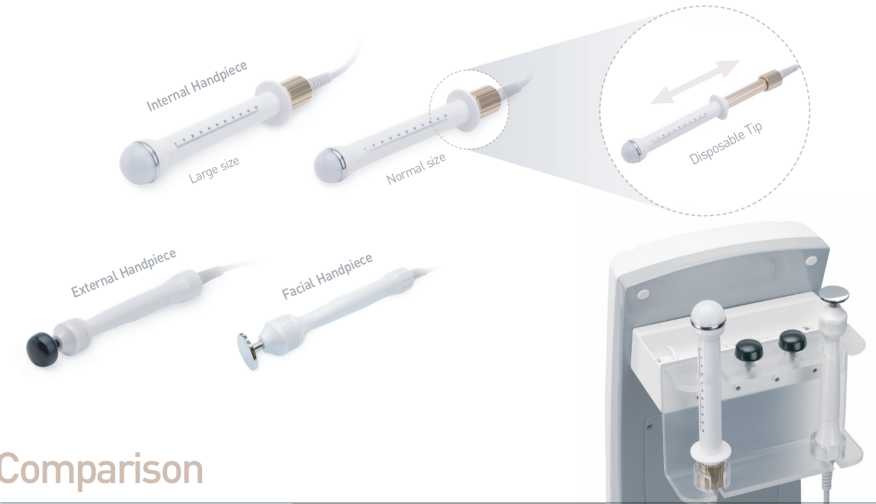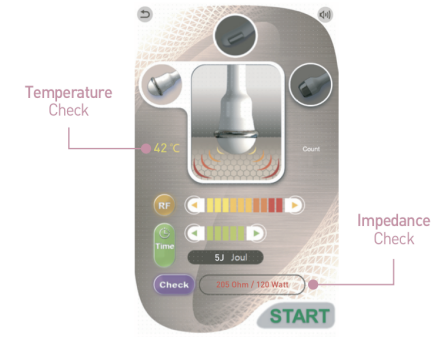

Comparison

	RF Energy Intima RF	Laser Procedure
Sexual Activity	The day of treatment	4-5 days after
External Treatment (Labia Remodeling)	Able	Unable
Smoke	None	Occurs to the point of disrupting the procedure
Bleeding	None	Occasional bleeding
Lubricant	No need	Need lubricant when inserting the handpiece
Moisture	Increase	Decrease
Recovery Time	No extra time needed	Around a week